


Getting to know your

Cochlear 'Nucleus 6' CP910 Processor

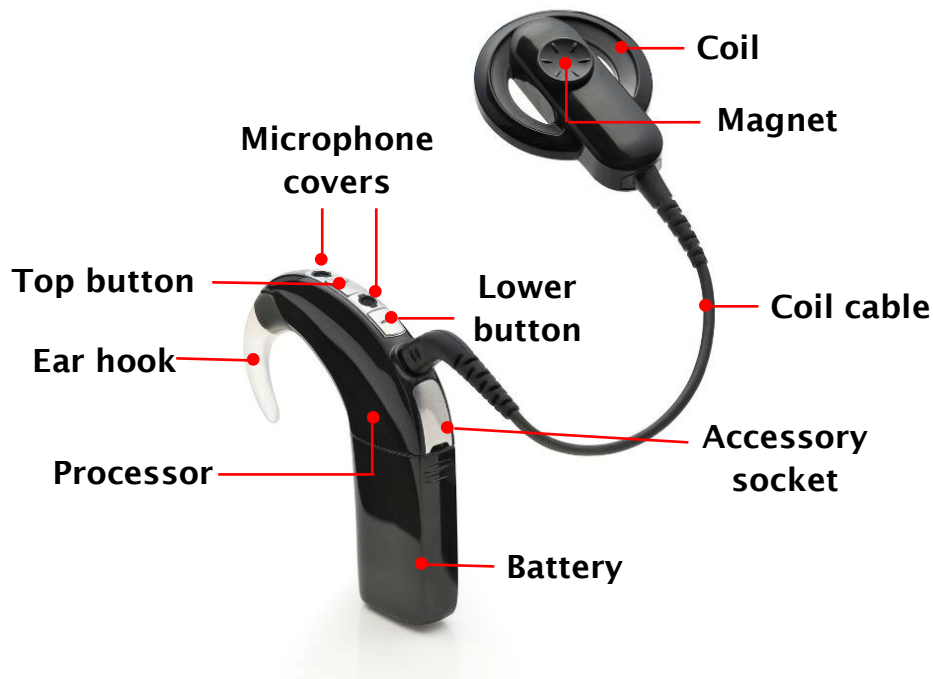

Switch 'on' by connecting the battery, 'off' by disconnecting the battery

TOP TIPS


- Change the microphone covers every 3 months
- To stop moisture affecting your processor, use your dry box every night
- If the processor sounds muffled, 'crackly' or stops working, change the cable. If the problem continues, contact the repairs service.
- If the processor looks damaged, contact the repairs service.
- The cochlear website can help you to troubleshoot any problems:

<http://www.cochlear.com/wps/wcm/connect/uk/home/support/cochlear-implant-systems/nucleus-6-system/troubleshooting/>

- Or visit the AIS website for helpful information and troubleshooting advice.

Getting to know your

CR230 remote control

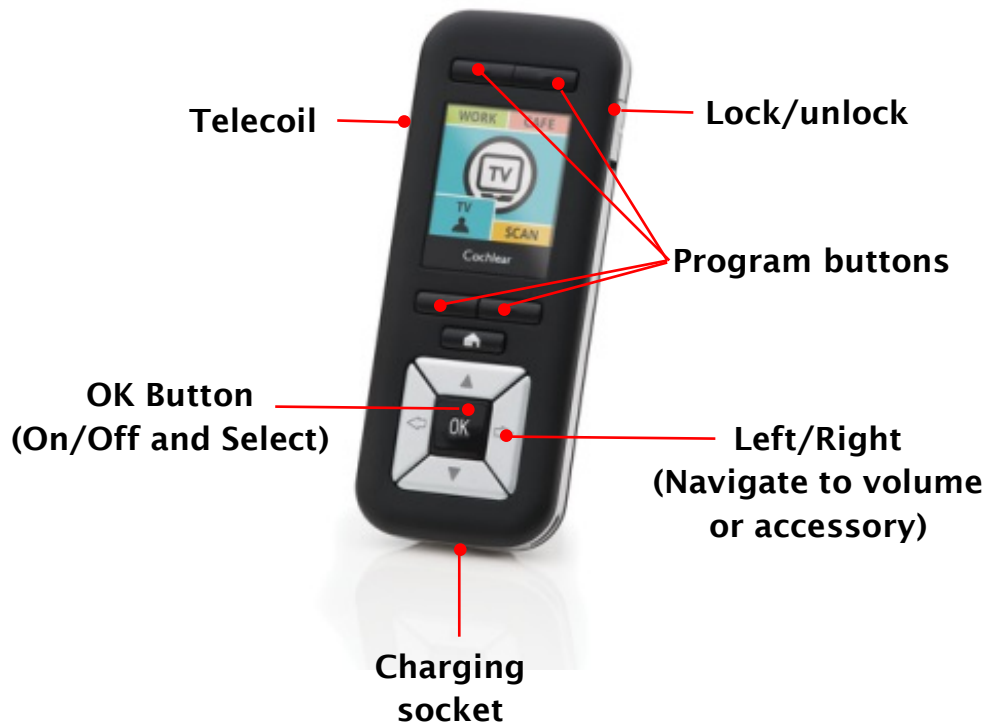

Switch your remote control 'on' or 'off' by pressing and holding the Ok button

TOP TIPS

- If you have a new processor you will need to pair it with the remote – to do this, place the processor coil on the back of the remote and press 'Ok' when prompted. To pair two processors to one remote, repeat the above step with your second processor.
- Check if your processor is working by quickly pressing the 'Ok' button on your remote – if the processor is on your head and working, the picture should be green in colour.
- You will need to charge your remote every few weeks using the mini USB cable in your kit.