


Getting to know your

Cochlear 'Nucleus 8' Processor

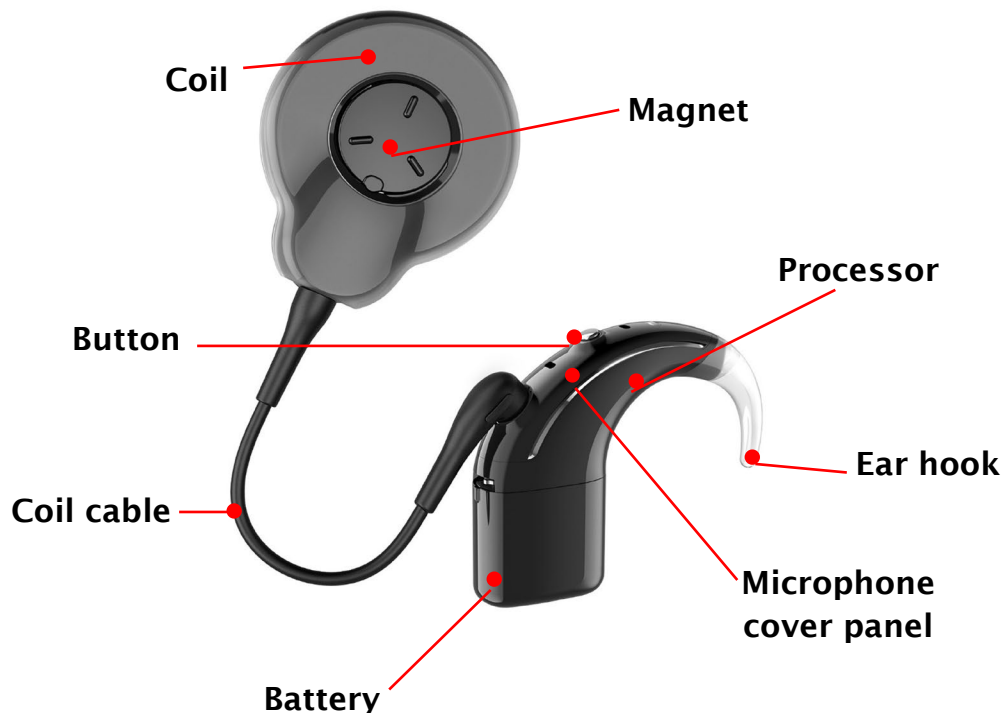

Switch 'on' by connecting the battery, 'off' by disconnecting the battery

TOP TIPS

- Change the microphone cover panel every 3 months.
- To stop moisture affecting your processor, use your dry box every night
- If the processor sounds muffled, 'crackly' or stops working, change the coil (you should have a spare in your kit). If the problem continues, contact Cochlear Care (or USAIS if you use us for repairs).
- If the processor looks damaged, contact the repairs service.
- The Cochlear website can help you to troubleshoot any problems:
<https://www.cochlear.com/uk/en/home/ongoing-care-and-support/device-support/cochlear-nucleus-portfolio/nucleus-8>
- Or visit the AIS website for helpful information and troubleshooting advice.

Getting to know your

CR210 remote control (not routinely offered with Nucleus 8)



 Switch your remote control 'on' or 'off' by sliding the side button.

TOP TIPS

- If you have a new processor you will need to pair it with the remote – to do this, place the processor coil on the back of the remote. To pair two processors to one remote, repeat the above step with your second processor.
- Your remote control does not need charging – it takes a large flat CR2032 battery that can be replaced by removing the back cover with the yellow screwdriver in your kit.
- The remote control automatically goes in to 'standby mode' to stop the battery draining. To 'wake it up' just switch it off and on again using the slider on the side.